




**27** CLINICAL EVALUATORS  
**660** TOTAL USES  
**86%** CLINICAL RATING

**Key features:** Universal bonding agent (light-cured) • Compatible with all etching modes • Multiple Indications

## Description

*HugeBond Universal* is a light-cured, single-component universal bonding agent which is:

- Used in self-etch, total-etch and selective-etch procedures
- Compatible with light-cured, dual-cured, and self-cured cements
- Reliable for bonding direct and indirect restorations
- Formulated with MDP monomer that enhances chemical adhesion to direct and indirect restorations such as dentin, enamel, composite, metal, and zirconia

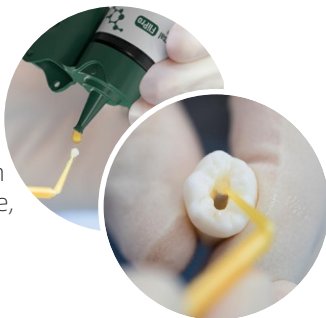


## Indications

- Direct restorative procedures
- Indirect restorative procedures
- As a primer for indirect restorations
- Repair of ceramic, composites, or metal restorations
- Desensitizer

DENTAL ADVISOR Biomaterials Test Results:

Shear Bond St rength Self-Etch Mode

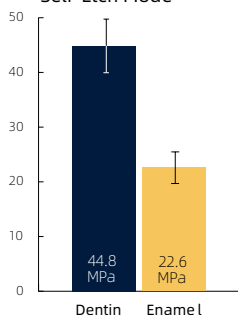

## Clinical Tip

- Use a small amount of air when thinning out the first coat



EASY TO USE WITH MANY INDICATIONS.

## Evaluators' Comments

"The ability to use a single bonding material for all styles and types of composites is very nice."

"The viscosity was thicker than most agents. I felt I had great coverage."

"The bottle design made it easy to get the amount you needed. Some bottle designs dispense more than what you need or want."

"Great consistency and clinical usability"

"Film thickness is very high."

"Material is light sensitive - be sure to cover before and during use."


88% Consultants who would: Recommend to a colleague

4% Consultants who would want to stock in office: Yes, instead of current product

62% Yes, in addition to current product

19% I might want to order this product for certain cases

### Evaluation Summary:

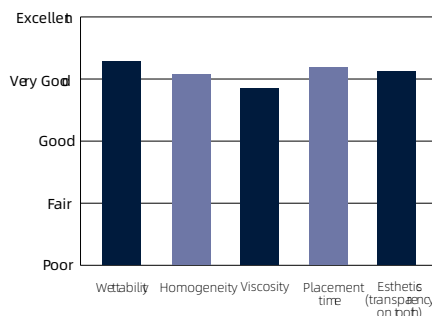

### Compared to Competitive Products:

